

Partecipating in children life





Partecipating in children life

TILT IN SPACE PAEDIATRIC FRAME

- SAFE, RESISTANT AND LIGHTWEIGHT: made in alloy of satin aluminium without welds
- HIGHLY MODULAR: able to follow the pathological evolution and the children growth
- VERSATILE: it can be used indoor and outdoor
- · EASY TO HANDLE
- · COLOURFUL





CADJUSTMENTS:


- CONTINUOUS DISPLACEMENT OF THE UPPER FRAME (which houses the positioning system) with respect to the underframe
- · SEAT DEPTH ADJUSTMENT
- HEIGHT FROM SEAT RAILS TO FLOOR (HIGH/LOW) ADJUSTMENT




- KNEE ANGLE
- DETACHABLE 75°/90°
 and ELEVATING






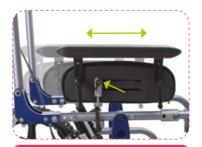

HORIZONTAL DISPLACEMENT

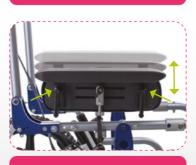

HEIGHT ADJUSTMENT

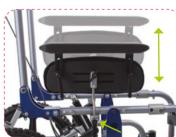

KIT EXTRACTION



WIDTH ADJUSTMENT

FUNDAMENTAL PARTS

- Frame with coloured plates made in alloy of satin aluminium
- · Adjustable single handlebar with folding backrest tubes
- Reclining mechanism by predetermined plate w/pin (3°,10°,18°)
- · Parking Brake lever
- Continuous tilt in space from 2° up to 35° (from 6° to 39°)

① Tilt in space can be up-graded of 4° with an increment of the starting inclination of the frame by operating on the positioning holes of the casters forks

- · Knee angle 90°
- Footplates
- Simple Calf Strap
- · Fixed anti-tip system





(WITH OR WITHOUT UP-CHARGE, PLEASE REFER TO THE ORDER FORM)

- · Casters:
- 150 mm or 6" poly
- 175 mm or 7" poly
- · Posterior wheels:
- 300 mm or 12" (poly or pneumatic, poly or pneumatic with drum brake, poly or pneumatic with drum brake + lever)
- 400 mm or 16" (poly, poly with drumbrake)
- 500 mm or 20" (poly or pneumatic without handrim, pneumatic with handrim and drum brake, pneumatic with handrim and drum brake + lever)
- · Knee angle detachable 75°, elevating
- · Footrest:
- Footboard
- Slotted only on footplates
- Tilt in space by pedal switch alternatively to the lever standard one, which in strongly recommended
- · Back posts:
- Single handlebar kit, reclining angle adjustable
- Telescopic push handles kit with lever block
- Fixed push handles kit
- · Quick Release pin for easy reclining
- · Detachable anti tip system
- Canopy
- · Foot muff
- · Rain Cover
- Basket

Note: for the complete list of the optional parts and/or accessories, please refer to the latest order form in force









HEIGHT (floor-to-seat-rails)

HEIGHT:

· Adjustable from 43 cm or 49 cm (400 mm wheels only 45 cm height)



ADACTA GIÒ IS AVAILABLE IN TITANIUM SATIN ANODISED WITH **3 DIFFERENT PLATES COLOURS**

BLUE	PINK	BORDEAUX				

FRAME SIZE

Frame Sizes

(outside tubes width x depth cm)

\	SEAT WIDTH	DEPTH AND ADJUSTMENTS
1	25	·32 ·34.5 ·37
2	25	• 36 • 38.5 • 41
3	30	·32 ·34.5 ·37
4	30	• 36 • 38.5 • 41
5	35	· 36 · 38.5 · 41
6	35	• 40 • 42.5 • 45
7	38	· 40 · 42.5 · 45











ADJUSTABLE HANDLEBAR KIT

MECHANICAL RECLINING

COLOURED PLATES

ANTI TIP SYSTEM









DEPTH ADJUSTMENT

SOLID BASE

DETACHABLE KNEE ANGLE

ELEVATING KNEE ANGLE

FRAME DIMENSIONS AND FEATURES

The device can be disassembled without the use of specific tools. It can be transported in the car or in the plane. It can be used as seat in a vehicle, provided that this one is equipped of transportation kit adaptations for people with disabilities.

LOAD (KG) / WEIGHT (KG)	1	2	3	4	5	6	7
	30/12.6	30/12.8	50/13.4	50/13.6	50/14.1	50/14.5	50/15.1



FRAME DIMENSIONS AND FEATURES

OVERALL DIMENSIONS (CM) OF THE FRAME WITH BACKREST TUBES (FIXED) FOLDED DOWN, FOOTREST AND POTERIOR WHEELS DETACHED									
SIZE	1	2	3	4	5	6	7		
LENGTH	64	66	64	66	66	68	68		
WIDTH	38	38	42	42	48	48	51		
HEIGHT	64	64	65	65	66	66	67		
WEIGHT(KG)	9.4	9.6	10.1	10.3	10.7	11.1	11.5		

Note: The length dimensions, referred to the average values, are subjected to variations due to the adjustments. The values reported are intended for the standard configuration.

WIDTH (CM)		1	2	3	4	5	6	7
300 MM WHEELS	standard	43	43	48	48	53	53	56
	w/drum brake	44	44	49	49	54	54	57
400/500 MM WHEELS	standard	45	45	50	50	55	55	58
	w/drum brake	48	48	53	53	58	58	61
PIVOT WIDTH (CM)		66	66	69	69	72	72	75
Q*(MIN-MAX) HEIGHT SEAT RAIL TO FOOTPLATE (CM)		13-20	13-20	13-23	13-23	13-23	13-26	13-26
H** SEAT RAIL TO FLOOR (CM)			43 and 49 with 300/500mm posterior wheels and 150/175mm castors 45 with 400 mm posterior wheel and 150/175mm castors					
MAX SL	MAX SLOPE 7°							

MAXIMUM LENGTH (CM) WITH 75° HANGER	1	2	3	4	5	6	7
300 MM WHEELS	86	89	86	89	89	92	92
400 MM WHEELS	87	90	87	90	90	94	94
500 MM WHEELS	91	94	91	94	94	97	97



Manufactured and Distibuted by:



Via Montagna Z.I. - Lotto 41
72023 Mesagne (BR) - ITALY
Tel. +(39) 0831 777840
Fax. +(39) 0831 730739
Email: sales@promedicare.it
www.versainserto.com
www.promedicare.eu

CERTIFIED COMPANY WITH QUALITY
MANAGEMENT SYSTEM ISO 9001 / EN ISO 13485







The producer reserves the right to make changes without notice. The contents, the form, the design, the image, the structure of this catalogue are proparty of Pro Medicare s.n.! Any improper or anauthorised use is strictly forbidden. Any improper or anauthorised use is strictly forbidden. Any imitation and/or explicit or not explicity copy can be pursued by the Law.